Vario Photometer
DP 300
Operating Manual
Version 5.13 / 5.13 SI
Edition 2021-01







Dear customer,

We are pleased that you have chosen the Vario Photometer from Diaglobal GmbH and thank you for the confidence you have placed in us.

The Vario Photometer belongs to a new generation of small mobile devices developed by Diaglobal GmbH and specially designed for on-site analysis.

With the software version V5.3 and higher, an automatic test of the device function has also been integrated. Therefore, the Vario Photometer complies with the requirements of the guidelines of the German Medical Association.

With the Vario Photometer, 12 clinical-chemical parameters can be determined. The device can be supplied with SI units of measurement on request (see chapter 9, Technical Data, table Measuring Ranges).

The kits and accessories required for the test are also available from Diaglobal GmbH.

All the best for your work with the new Vario Photometer!

Yours Diaglobal GmbH

Table of contents

		Page
1.	General information regarding the Photometer	4
2.	Installation	5
3.	Description of the device	5
3.1	Power supply	6
3.1.1	Mains power operation	6
3.1.2 3.2	Mains-independent operation Measuring system	6 6
4.	Service	7
4.1	Adjustment and Calibration	7
4.2	Maintenance	7
4.3	Cleaning Instructions	7
4.4	Malfunctions	7
4.5	Disposal	7
5.	Required reagents and laboratory accessories	8
5.1	Expiration date of consumables	8
5.2	Reagents / parameter list	8
5.3	Control materials	9
5.4	Laboratory aids and accessories	9
6.	Quality control according to the Guideline of the German Medical Association	10
7.	Measuring procedure	11
7.1	End-point measurement	11
7.2	End-point measurement, taking the sample's blank value and	
	pre-programmed measuring time into account	11
7.3	Multipoint measurement, taking the sample's blank value and	12
7.4	recognition of the end point into account Multipoint measurement, taking the sample's blank value and	12
7.4	calculation of the end point into account	12
8.	Measurement	13
8.1	Switching the device on	13
8.2	Self-test when switching on	13
8.3	Test selection	13
8.4	Switching the device off	13
8.5 8.6	Integrated operational device checks	13 14
0.0	Notes on taking samples and carrying out measurements	14
9.	Technical data	16
10.	General guidelines and notes	17
11.	Appendix: "Step-by-step measurement"	17ff.

1. General information on the Photometer

Device name: Vario Photometer

Model: DP 300

Features: In-vitro diagnostics, measuring device for the

determination of selected clinical-chemical parameters in blood, serum/plasma and

cerebrospinal fluid

The Vario Photometer fulfils the basic requirements of Appendix I of Directive 98/79/EC regarding in-vitro diagnostics.

The conformity of the device with Directive 98/79/EC is confirmed by the use of the CE marking.

Manufacturer: Diaglobal GmbH

Innovationspark Wuhlheide Köpenicker Str. 325 / Haus 41

12555 Berlin

Tel: +49 (0) 30 6576 2597 Fax: +49 (0) 30 6576 2517 E-Mail: info@diaglobal.de

http://www.diaglobal.de

2. Installation

For trouble-free operation of the device, the following environmental conditions must be met:

- Ambient temperature: 0 °C ... 40 °C
- No direct exposure to sunlight or similar sources of radiant heat
- Free from excessive dust
- Free from vibrations
- Free from interference by electromagnetic waves
- Operation on a horizontal surface

Please observe the following instructions for use:

Insert a rechargeable battery or normal battery if the device is to be operated independently of a power supply or connect the photometer to a power supply unit.

Press the **<ON/ENTER>** key (Fig. 1) to activate the internal device check which is automatically carried out by the device.

The device is then immediately ready for measurement.

3. Description of the device



Fig. 1

3.1 Power supply

The Vario Photometer can be operated as desired using a power supply, a (9V block) battery or (model 6F22 or PP3) rechargeable battery.

3.1.1 Mains power operation

The Photometer is supplied with a power supply unit for operation on a mains voltage in the range of 100 V ... 240 V AC. The mains plug is marked with a Diaglobal logo (sticker).

The connector plug of the power supply unit is connected to the power supply socket on the back of the device.

3.1.2 Network-independent operation

To insert the rechargeable battery or the normal battery:

Unscrew the knurled screws on the bottom of the unit and remove the battery compartment cover. Connect the battery to the push-button contact and insert it into the device. Replace the battery compartment cover and screw in the knurled screws.

Please note:

The Vario Photometer can be operated using a power supply without the need to remove the rechargeable battery or the normal battery.

The rechargeable battery cannot be charged while it is installed. A separate battery charger is required for this purpose.

3.2 Measuring system

The optical section is shown in Fig. 2.

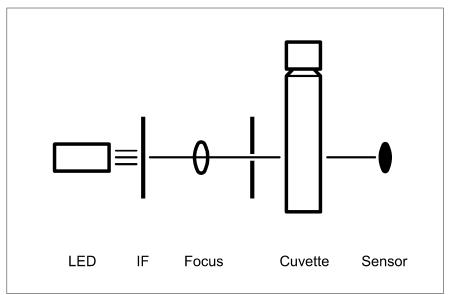


Fig. 2

The light emitted by an LED is first selected into its wavelength ranges (520 nm and 546 nm) by an interference filter IF (HBW \sim 5 nm) and then bundled and directed onto the cuvette in the shaft. After passing through the cuvette, a broadband photosensor converts the light falling on its sensor surface into a current, proportional to the intensity.

4. Service

4.1 Adjustment and Calibration

The instrument is adjusted and calibrated at the factory on delivery, adjustment by the customer is not necessary.

Adjustment is carried out via the interface socket on the rear panel. It can only be carried out at the factory, adjustments by the customer are not possible.

Information on calibrating the device can be found in chapter 6. *Quality control according to the Guideline of the German Medical Association*.

4.2 Maintenance

The device is maintenance-free. Maintenance after the warranty period is recommended, but not mandatory.

Due to the integrated test of the unit functions (chapter 8.5) and regular tests with control material, maintenance is only recommended if one of these two test functions indicates an error message.

4.3 Cleaning Instructions

Commercially available decontaminating solutions commonly used in clinical chemistry laboratories, such as Mikrozid® AF Liquid, Bacillol® plus, 3 % Kohrsolin® or similar, are recommended for cleaning the device and the surface. Before cleaning the unit with a soft cloth and the decontaminating solution, it must be switched off and the electrical power supply must be disconnected.

Make sure that no liquids get into the device. There is no protection against penetrating liquids (Code IP X0).

The cuvette shaft must not be cleaned by the user of the device, as this may damage the device. If cleaning is necessary, especially because of leaking liquids or broken glass, please contact Diaglobal GmbH.

4.4 Malfunctions

If any malfunctions or problems occur, simply call us. Most questions can be answered on the phone. Non-functional units should be sent to our Berlin address. We will provide a loan device for the duration of the repair.

4.5 Disposal

Diaglobal GmbH will take back and dispose of units that are no longer needed or cannot be repaired, free of charge.

5. Required reagents and laboratory accessories

5.1 **Expiration date of consumables**

It is important to ensure that all consumables may only be used within the expiration date.

5.2 Reagents / parameter list

The following tests can be measured with the Vario Photometer:

Doubles	Sample material			Toots (no sk	Aut no	
Parameter	Blood	Serum	Plasma	Tests/pack	Art. no.	
Lactate	+	-	+	40	LAC 142	
Lactate-rapid	+	-	+	40	LAC 342	
Glucose	+	+	+	40	GLU 142	
Cholesterol	+	+	+	40	CHO 142	
HDL-Cholesterol ^{1) 3)}	+	+	+	20	HDL 321	
Triglycerides	+	+	+	40	TRI 142	
Haemoglobin (HiCN-method)	+	ı	-	40	HB 142	
Haemoglobin (SLS-method)	+	-	-	40	HB 342	
Erythrocytes	+	-	-	40	ERY 142	
Haematocrit	+	-	-	40	HCT 142	
Urea ^{1) 3)}	+	+	+	20	HST 321	
Uric acid	-	+	+	40	HSR 342	
Creatinine ²⁾	-	+	+	20	KRE 121	
Bilirubin	-	+	+	40	BIL 142	
Neonatal Bilirubin ^{1) 3)}	+	+	+	40	DIL 142	

Mini centrifuge required (Art. no. DZ 002)
 Dry Block Thermostat required (Art. no. DZ 003)
 From blood, with subsequent sample preparation (centrifugation with mini centrifuge)

5.3 Control materials

Art. no.	Description	Contents
HEM QS	Haemoglobin control Haemolysate for correctness and precision control of haemoglobin determination in blood in the normal range	5 x 1 mL
ERY QS	Erythrocytes- and Haematocrit control Control blood for accuracy and precision control of haematocrit and erythrocytes determination in blood in the normal range	5 x 1 mL
GLU QS	Glucose control 100 mg/dL	3 x 4 mL
LAC QS	Lactate control set 2 mmol/L; 4 mmol/L; 10 mmol/L	3 x 4 mL
BIL QS	Bilirubin control caps Lyophilizates for accuracy and precision control of the bilirubin measurement	20 caps

5.4 Laboratory aids and accessories

Art. no.	Description	Contents
LH 001	Blood lancets	500
LH 004	Capillaries 10 μL, with ring mark	250
LH 006	Cuvette rack	1
LH 007	Micropipettor (pipetting aid)	1
LH 009	Cellulose swabs	500
LH 010	Cellulose swab box	1
LH 011	Alcohol pads, non-sterile	100
LH 012	Powder-free nitrile gloves size M	200
LH 020	Capillaries 20 μL, heparinised, end-to-end	100
LH 023	Capillaries 100 μL, with ring mark	250
LH 024	Capillaries 20 µL, with ring mark	250
LH 031	Microvettes CB 300	100
LH 034	Safety Lancets Neonatal, pink 1.2 mm	200
LH 035	Safety lancets extra, orange 1.8 mm	200
LH 050	Reaction tubes to separate the plasma	500
LH 055	Pipette tips 50-1000 μL blue, for pipette LH 500	500
LH 060	Capillaries 60 μL, heparinised, end-to-end	5 x 20
LH 500	Pipette fix 500 μL	1
DZ 002	Mini centrifuge	1
DZ 003	Dry Block Thermostat	1

All reagent kits, control materials and other materials are supplied by Diaglobal GmbH and can be stored and transported together with the Vario Photometer in a practical case.

6. Quality control according to the Guideline of the German Medical Association¹⁾

The Vario Photometer has been specially developed for near-patient immediate diagnostics with unit-use reagents (German Medical Association, part B, chapter 2.1.5). According to the guideline of the German Medical Association, there is therefore no obligation to participate in surveys (German Medical Association, part B, chapter 2.2, paragraph (3) a). The user only has to carry out internal quality checks.

Internal quality assurance is carried out in the form of a weekly accuracy check (calibration) with subsequent documentation of the measured value. The corresponding protocol forms are available from Diaglobal free of charge. We recommend using the Diaglobal control solutions LAC QS and GLU QS to check the accuracy of lactate and glucose determinations.

We recommend using the blood control HEM QS and ERY QS with target values in the normal concentration area for checking the accuracy of determinations of haemoglobin, haematocrit and erythrocyte counts.

We recommend using the Diaglobal control caps BIL QS to check the accuracy of bilirubin determinations.

For all other parameters we recommend using the universal control sera from the company, Roche, www.roche.de:

PreciControl ClinChem Multi 1 Order-No.: 05 947 626 190 (4 x 5 mL) for normal range PreciControl ClinChem Multi 2 Order-No.: 05 947 774 190 (4 x 5 mL) for pathological range

In agreement with the requirements of the German Medical Association, a test of the device function (see operating instructions, chapter 8.5) is integrated in the Vario Photometer, therefor a daily test by means of a standard manual test (German Medical Association, part B, chapter 2.1.5, paragraph (2) is not necessary.

The Vario Photometer is suitable for the speedy detection of gestational diabetes and fulfils the requirements of the Maternity Guidelines²⁾ and the S3-Directive³⁾. Glucose can be measured from whole blood as well as from venous plasma. The displayed measured value is - according to the requirements - always related to venous plasma.

¹⁾ Guideline of the German Medical Association for the quality assurance of laboratory medical examinations Deutsches Ärzteblatt | Jg. 116 | Heft 51-52 | 23. Dezember 2019

²⁾ BAnz. Nr. 36, S914

³⁾ AWMF-Register Nr. 057/008

7. Measuring process

7.1 Endpoint measurement

The absorbance is measured after reaching the endpoint.

It is measured against the reagent's blank count.

Parameters: Parameters: Haemoglobin SLS (HB SLS), Erythrocytes (ERY), Haematocrit (HCT), Haemoglobin (HiCN), Uric Acid (UA), HDL Cholesterol (HDL/CHO)

Calculation: Concentration = Absorbance x Factor

The erythrocyte and haematocrit counts are determined using stored reference curves.

7.2 Endpoint measurement with consideration of the sample blank value and pre-programmed measuring time

After measuring the sample blank value, the colour reaction in the cuvette is started and the endpoint absorbance is measured after a specified time has elapsed.

Parameter: Creatinine (CRE)

Calculation: Concentration = Absorbance x Factor

Measuring time: 2 minutes

The samples are measured one after the other: Sample 01: Measurement 1 (sample blank value)

Sample 01: Measurement 2 (result)

Sample 02: Measurement 1 (sample blank value)

Sample 02: Measurement 2 (result)

etc.

Parameter: UREA

Calculation: Concentration = Absorbance x Factor

Measurement time: 10 minutes

Blank values of the samples are measured one after the other:

Sample 01: Measurement 1 (sample blank value) Sample 02: Measurement 1 (sample blank value) Sample 03: Measurement 1 (sample blank value)

etc.

Results are measured one after the other:

Sample 01: Measurement 2 (result) Sample 02: Measurement 2 (result) Sample 03: Measurement 2 (result)

etc.

7.3 Multi-point measurement with consideration of the sample blank value and recognition of the endpoint

After measuring the sample blank value (=measurement 1) the colour reaction in the cuvette is started. The reaction process is monitored by the device (=measurement 2). The measuring process is terminated as soon as the endpoint is reached.

The time needed to reach the endpoint depends on the temperature. It is normally 2 - 6 minutes for the lactate test. If temperatures are close to freezing point, measuring times can take up to 20 minutes, depending on the parameters. You can choose between single and series measurements up to a maximum of 20 samples.

For single measurements, the samples are processed one after the other. For series measurements, all A1 values are measured first.

Parameters: Lactate (LAC), Cholesterol (CHO), Triglycerides (TRI) Calculation: concentration in plasma = Absorbance Difference x Factor

7.4 Multi-point measurement with consideration of the sample blank value and calculation of the endpoint

After measuring the sample blank value (=measurement 1) the colour reaction in the cuvette is started. The course of the reaction is monitored by the instrument. The endpoint is calculated using several absorbance values recorded at different times.

Parameters: Glucose (GLU), Lactate-rapid (LAC-rapid)

Measuring times: Glucose 2 minutes Lactate-rapid 1 minute

8. Measurement

8.1 Switching the device on

Press the **<ON/ENTER>** key

8.2 Self-test when switching on

When the device is switched on, a self-test of the digital and analogue circuitry is conducted. The operational device check proceeds automatically after it is switched on. It takes approx. 5 seconds, after which the unit is ready for measuring.

Note:

If it becomes obvious during the test that one of the device functions does not correspond to the required settings, <SERVICE> will appear in the display.

In this case, switch the device off.

Please call Diaglobal GmbH service (Tel. +49 (0) 30 6576 2597) or contact your specialist retailer.

8.3 Test selection

Press the **<ON/ENTER>** key.

The desired test is selected from the menu with the right or left arrow key:

```
LAC - LAC-rapid - GLU - CHO - HDL/CHO - TRI - HB - HB-SLS - ERY - HCT - UREA - UA - CRE - BIL - BIL N - ABS520 - ABS546
```

Pressing the right arrow key activates the next test while pressing the left arrow key returns to the previous test. The selected test is shown in the upper right corner of the display.

Confirm test selection with the **<ON/ENTER>** key.

8.4 Switching the device off

To switch the device off, press both arrow keys simultaneously.

8.5 Integrated operational device checks

Self-test when switching on

Testing of the digital and analogue circuits of the device is automatically performed by the device when it is switched on. Please see chapter 8, point 8.2.

<u>Differential measurements</u>

All measurements are based on differential measurements. I. e. after selecting the desired test, the device requests a zero measurement with a blank value cuvette. This creates a reference base to the measured value so that minor deviations can be compensated.

Measurement range controls

The measurement ranges of all measurement results shown in the display are verified by an integrated measurement range control. If the measurement range is exceeded, an error is displayed.

The measurement ranges that are separately defined for each parameter are documented on the respective package inserts as well as in this operating manual, chapter 9, Technical Data.

Plausibility controls

For multi-point measurements, the absorbance measured first forms the reference basis. The programme verifies the plausibility of the individual measured values. If specific requirements (e.g. E2 > E1 for ascending reactions) are not met, an error message is displayed.

8.6 Notes on taking samples and carrying out measurements

Errors in taking samples will always lead to incorrect measurement results. This chapter addresses the most common errors that can occur during taking samples and measuring samples.

- 1. Before measuring, cuvettes stored in a refrigerator must be brought to room temperature. If the cuvettes are too cold, they will become misty with water on the outer wall due to the humidity, which will lead to incorrect measurement results.
- 2. Never touch the lower part of the cuvette (where the liquid is) with bare hands. If this should happen accidentally, clean the vials with a fluff-free cloth before use. Cleaning with a fluff-free cloth is recommended in any case. Even if the package is still new and unopened. Fingerprints on the cuvette lead to incorrect measurement results.
- 3. If the blood is taken from the heel using the microvette (neonatal bilirubin), make sure that there is enough blood (approx. 60 μ L) in the microvette, as 20 μ L serum/plasma is needed for the measurement. Close the microvette well after taking the blood sample and return it to the small sample tube for centrifugation.
- 4. After centrifugation of the microvette, please make sure that the centrifugate has separated completely and that the supernatant is clear and free of any solid particles. If not, repeat the centrifugation. If the supernatant is not free of suspended matter or if particles of the centrifugate accidentally enter the capillary, the measurement result will be incorrect.
- 5. If blood is taken from the fingertip or earlobe, note that the first drop that forms spontaneously must be wiped away with a cellulose swab. It contains a high proportion of tissue fluid, which will corrupt the measurement result.
- 6. The second drop that forms is for blood sampling. To support blood collection, it may be pressed carefully (!). The emphasis on carefully, otherwise too much tissue fluid will get into the blood sample again.
- 7. Make sure that the blood drop that forms is large enough to fill the capillary with the required sample volume in one go. Repeated filling of the capillary leads to air bubbles that cannot be removed from the capillary. If air bubbles form, discard the capillary and start sampling again.

8. The capillary must be filled exactly up to the black ring mark.

Please note: A deviation of only 1 mm from the ring mark is sufficient to obtain a completely incorrect measurement result!

If the sample is above the black ring mark, this will lead to incorrect positive measurement results. A cellulose swab can be used to carefully soak up too much blood.

If the sample is below the black ring mark, this will lead to incorrect negative measurement results. In this case, correction is hardly possible due to the air bubble that will form when trying to collect more blood.

- 9. Before the capillary is placed in the cuvette, the lower area must be carefully wiped on the outside with a cellulose swab to remove sample particles attached to the capillary. Otherwise, this would lead to incorrect positive measurement results.
- 10. With the help of the micropipetter, the sample is completely transferred into the cuvette. The complete transfer of the sample is done by ejecting it several times with the help of the push button on the micropipetter.

Please note: The micropipetter is only used when the capillary is filled with the sample. It is not needed for filling the capillary. The capillary is filled by the capillary action alone.

- 11. During series measurements, make sure that the order of the samples is not reversed. Otherwise, the device cannot assign the samples correctly, which leads to unreliable measurement results.
- 12. When changing the cap with the starter cap, make sure that the substance in the starter cap has completely dissolved. Failure to do so will result in a non-linear kinetic reaction process, which will lead to an error message in the display or unreliable measurement results.

9. Technical data

Storage temperature: $-20 \, ^{\circ}\text{C} \dots 70 \, ^{\circ}\text{C}$ Operating temperature: $0 \, ^{\circ}\text{C} \dots 40 \, ^{\circ}\text{C}$

Dimensions: $200 \times 100 \times 50 \text{ mm}$

Weight: 450 g

Measuring principle: Absorption measurement with single beam

photometer (Fig. 2), chopped operation

Projector: LED

Spectral apparatus: Interference filter

Measuring wavelengths: 520 nm and 546 nm

Spectral half-width value: $\sim 5 \text{ nm}$ External light influence: Negligible

Interface: V24 (9600, 8, n, 2) Power supply: $6 \text{ V} \dots 12 \text{ V} \text{ DC}$ Current consumption: max. 250 mA

Warm-up time: 0 min

Interference suppression: According to DIN VDE 0871 and DIN VDE 0875

Inaccuracy: < 0.5 % at A = 1.000

Relative photometric

short-time standard deviation: < 0.1 %

Measuring ranges:	DP 300	<u>DP 300 SI</u>	
Lactate	0.2 - 30 mmol/L	0.2 - 30 mmol/L	
Lactate-rapid	0.2 - 20 mmol/L	0.2 - 20 mmol/L	
Glucose	20 - 630 mg/dL	1.1 - 35 mmol/L	
Cholesterol	20 - 1300 mg/dL	0.5 - 35 mmol/L	
HDL-Cholesterol	10 - 200 mg/dL	0.2 - 5 mmol/L	
Triglycerides	20 - 2000 mg/dL	0.2 - 23 mmol/L	
Haemoglobin (HiCN-method)	0.0 - 50 g/dL	0.0 - 31 mmol/L	
Haemoglobin (SLS-method)	0.0 - 50 g/dL	0.0 - 31 mmol/L	
Erythrocytes	1.0 - 10 Mio/μL	1.0 - 10 Mio/μL	
Haematocrit	5 - 90 %	5 - 90 %	
Urea	5 - 200 mg/dL	0.8 - 35 mmol/L	
Uric acid	0.0 - 25 mg/dL	0.0 - 1500 μmol/L	
Creatinine	0.0 - 5 mg/dL	0.0 - 440 μmol/L	
Bilirubin	0.50 - 25 mg/dL	8.50 - 428 μmol/L	
Neonatal Bilirubin	2.30 - 50 mg/dL	39.0 - 850 μmol/L	
ABS 520 nm	A = 2.500	A = 2.500	
ABS 546 nm	A = 2.500	A = 2.500	

10. General Guidelines and Notes

EC Directives

1. Directive 98/79/EC on in-vitro diagnostic devices

EN / ISO standards

- **2.** EN ISO 9001:1994, Quality Management Systems, Model for quality assurance in design, development, production, installation and customer service
- **3.** EN ISO 13485, Medical devices, Requirements for regulatory purposes (application of EN ISO 9001)
- **4.** EN ISO 14971, Medical devices Application of risk management to medical devices
- **5.** EN 61010 -1, Safety requirements for electrical equipment for measurement, control and laboratory use Part 1: General requirements
- **6.** EN 61010 -2-101, Safety requirements for electrical equipment for measurement, control and laboratory use Part 2-101: Particular requirements for in-vitro diagnostic (IVD) medical equipment
- **7.** EN 61326 -1, Electrical equipment for measurement, control and laboratory use EMC requirements Part 1: General requirements
- **8.** EN 61326 -2-6, Electrical equipment for measurement, control and laboratory use EMC requirements Part 2-6: Particular requirements In-vitro diagnostic (IVD) medical equipment
- **9.** EN 592, Instructions for use for in-vitro diagnostic instruments for professional use

National directives and recommendations (Germany)

10. Guidelines for Quality Assurance of Laboratory Examinations of the German Medical Association of 23.12.2019

Note on electromagnetic compatibility

- a) The photometer meets the requirements for electromagnetic radiation and interference immunity as described in the IEC 61326 series of standards.
- b) Do not use this device near sources of intense electromagnetic radiation because they may interfere with correct functioning. A distance of at least 1 m should be maintained between an operational (switched on) mobile phone and the photometer during measurement.

Note on the unit's internal quality control

The functionality of the device is checked when it is switched on. In addition, electronically controlled checks are carried out for individual tests during the measurement, which leads to an error message if specified requirements are not met.

11. Appendix: "Step-by-step measurement"

Please refer to the illustrations in the "Step by step" instruction manual.

Device manual





1. Switch on: Press ON/ENTER key Wait for device check and confirm with ON/ENTER



2. Select test: Press arrow key until required test appears




3. Confirm required test: Press ON/ENTER



4. Switch off: Press both arrow keys at the same time

Note: If SERVICE appears in the display after the device check, the device has a defect. In this case, please contact our customer service at +49 (0) 30 6576 2597.

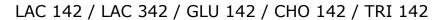



Single measurement



1. Insert capillary with 10 μL sample into cuvette



2. Eject sample several times with micropipetter into cuvette




3. Screw cap on Turn cuvette upside down several times



4. Switch photometer on with ON/ENTER key Wait for device check, confirm with ON/ENTER

Select the required test, confirm with ON/ENTER

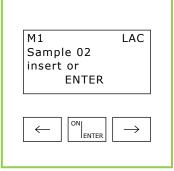

5. Zero point adjustment: Insert cuvette with sample (Fig. 3) into photometer, zero point is stored in memory


Remove cuvette after signal tone



6. Replace screw cap with starter cap



7. Turn cuvette upside down several times

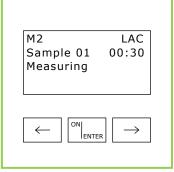

8. First press ON/ENTER Then insert cuvette into photometer



9. Time is displayed, wait for measured value

LAC 142 / CHO 142 / TRI 142




Number of samples per series: Up to 20 samples at the same time




1. Eject all samples one after the other several times with micropipetter into cuvette



2. Screw all caps on again Turn cuvettes upside down several times

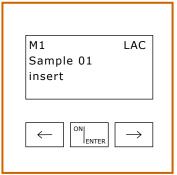

3. Switch photometer on with ON/ENTER key
Wait for device check, confirm with ON/ENTER
Select the required test, confirm with ON/ENTER

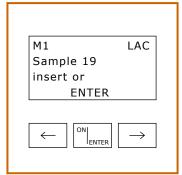

4. Zero point adjustment: Insert cuvettes with samples (Fig. 2) one after the other into photometer, all zero points are stored in memory

Note: Ensure the correct order


of the samples!



5. After the zero point adjustment of the last cuvette replace all screw caps with starter caps



6. Turn all cuvettes **simultaneously** upside down, repeat several times

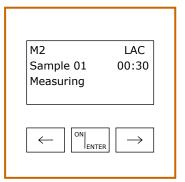

7. First press ON/ENTER key **Then** insert 1st cuvette into photometer


Time is displayed, wait for measured value



8. Read the measured value of the 1st cuvette, remove cuvette

Insert 2nd cuvette, read the measured value, remove cuvette,

and so on

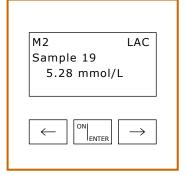

9. Insert the last cuvette, read the measured value, remove cuvette

Note: Ensure the correct order of the samples!





Number of samples per series: Up to 20 samples at the same time

Additionally required: Mini centrifuge, Cholesterol CHO 142



1. HDL 321 Transfer 60 μ L of the sample with an end-to-end capillary into each reaction tube "R" and mix well

Leave it for 5 minutes



2. HDL 321

Insert reaction tube "R" with capillary into centrifuge

Centrifugate for 5 minutes


Continue with CHO 142




Insert capillary with sample into cuvette



4. CHO 142
Eject sample several times with micropipetter into cuvette
Screw cap on
Turn cuvette upside down several times

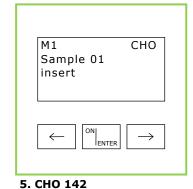

Switch photometer on with ON/ENTER key, wait for device check and confirm with ON/ENTER Select HDL/CHO, confirm with ON/ENTER Start with the measurement of CHO

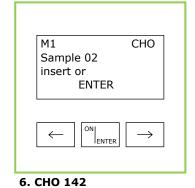

Zero point adjustment: Insert cuvette with sample (Fig. 4) into photometer, zero point is stored in memory


Remove cuvette after signal tone



7. CHO 142Replace screw cap with starter cap of CHO 142



8. CHO 142Turn cuvette upside down several times

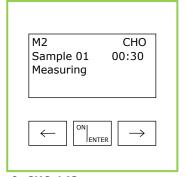

9. CHO 142
First press ON/ENTER
Then insert cuvette into photometer
Time is displayed, wait for measured value





Number of samples per series: Up to 20 samples at the same time

Additionally required: Mini centrifuge, Cholesterol CHO 142




10. CHO 142 The cholesterol value is stored in memory

Now continue with the measurement of HDL 321



11. HDL 321
Pipette 500 µL supernatant from the centrifuged reaction tube "R" (Fig. 2) into the cuvette


HDL 321

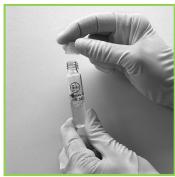


12. HDL 321Screw starter cap of HDL 321 on



13. HDL 321Turn cuvette upside down several times
Wait for 5 minutes

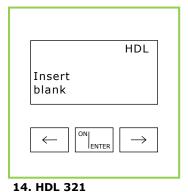

Zero point adjustment: Take an unprocessed HDL 321 cuvette (blank) from the package and insert it into the photometer Zero point is stored in memory

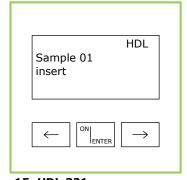

15. HDL 321 Remove blank after signal tone



16. HDL 321Insert cuvette with sample (Fig. 13) in photometer
Read measured value

Note: In order to carry out a series measurement, all CHO 142 values must be measured first

All CHO 142 values will be stored in memory, one after the other

Important: Pay attention to the correct order and assignment of the samples!








1. Insert capillary with 10 μ L blood sample into cuvette




2. Eject blood several times with micropipetter into cuvette




3. Screw cap on Turn cuvette upside down several times Wait 3 minutes HB 342: wait 30 seconds

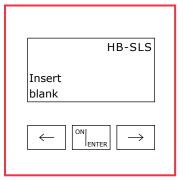


4. Switch photometer on with ON/ENTER key
Wait for device check, confirm with ON/ENTER
Select the required test, confirm with ON/ENTER

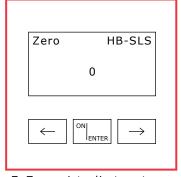

5. Zero point adjustment: Insert an unprocessed cuvette from package into photometer
Zero point is stored in memory

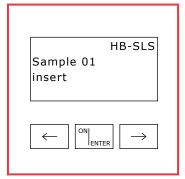

6. Remove cuvette after signal tone

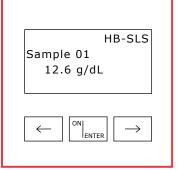

7. Insert cuvette with blood sample (Fig. 3) into photometer Read measured value



measurement:After zero point setting any number of additional samples can be measured

HST 321

diaglobal

Number of samples per series: Up to 20 samples at the same time

Additionally required: Mini centrifuge




1. Insert 20 µL sample with an end-to-end capillary into the reaction tube "R" and mix well




2. Insert reaction tube "R" with capillary into mini centrifuge
Centrifugate for 1 minute
Note: Ensure an even loading inside the mini centrifuge




3. Pipette 500 μL supernatant from the reaction tube "R" into the cuvette

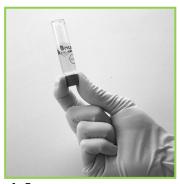

4. Screw cap on Turn cuvette upside down several times

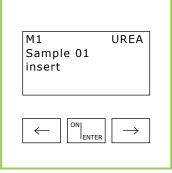

5. Switch photometer on with ON/ENTER key Wait for device check, confirm with ON/ENTER Select UREA, confirm with ON/ENTER

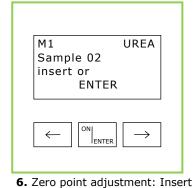


cuvette with sample (Fig. 4) into photometer, zero point is stored in memory

Remove cuvette after signal



7. Replace screw cap with starter cap Turn cuvette upside down several times

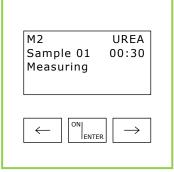

8. First press ON/ENTER Then insert cuvette into photometer Time is displayed

9. Wait for measured value (10 minutes)
Note: Series measureme

Note: Series measurement of Urea up to 20 samples is the same procedure as the serial measurement of LAC 142





1. Insert capillary with 100 μL sample (serum/plasma) into

Eject sample several times with micropipetter into cuvette




2. Replace screw cap with starter cap




3. Turn cuvette upside down several times




4. Leave the cuvette for 10 minutes at room temperature

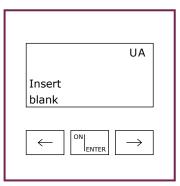

5. Switch photometer on with ON/ENTER key
Wait for device check, confirm with ON/ENTER
Select UA test, confirm with ON/ENTER

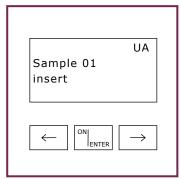

6. Zero point adjustment: Insert an unprocessed cuvette from package into photometer Zero point is stored in memory

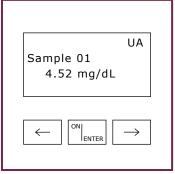

7. When the waiting time of 10 minutes has passed, insert the cuvette (Fig. 4) into the device Read measured value



In regard to series measurement:

After zero point setting any number of additional samples can be measured

KRE 121



Single measurement

Additionally required: Dry block thermostat (30 minutes preheated)



1. Pipette 1 mL buffer solution into each cuvette, then close cuvettes

Note: The example shows 6 samples. They can only be processed one after the other, serial measurement is not possible



2. Incubate cuvettes in the preheated dry block thermostat for 7 minutes




3. Transfer 500 µL of the sample (serum/plasma) into the cuvette

Then insert the cuvette **immediately** into the dry block thermostat




4. Incubate the cuvette with the sample for exactly 1 minute

Already during the incubation time, switch the photometer on with ON/ENTER

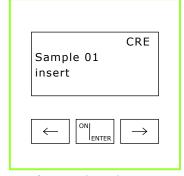

5. After switching the photometer on, wait for device check and confirm with ON/ENTER

Select CRE, confirm with ON/ENTER

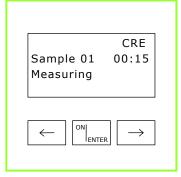

6. Insert cuvette into photometer

Time is displayed



7. Wait for measured value (2 minutes)



8. Read the measured value, remove the cuvette and repeat the same process with all other cuvettes, starting from Fig. 3

BIL / BIL N (BIL 142)




Additionally required: Mini centrifuge




1. Insert capillary with sample (serum/plasma) into cuvette BIL: 100 µL (adults) BIL N: 20 µL (newborns) Sampling of BIL N see next page




2. Eject sample several times with micropipetter into cuvette

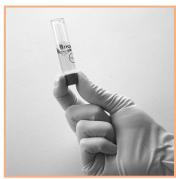


3. Screw cap on Turn cuvette upside down several times

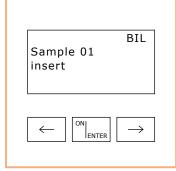

4. Switch photometer on with ON/ENTER key
Wait for device check, confirm with ON/ENTER
Select the required test, confirm with ON/ENTER

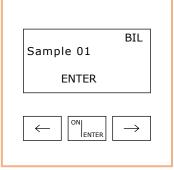


5. Zero point adjustment: Insert cuvette with sample (Fig. 3) into photometer, zero point is stored in memory Remove cuvette after signal tone

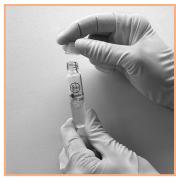

6. Replace screw cap with starter cap



7. Turn cuvette upside down several times

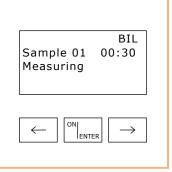

8. First press ON/ENTER Then insert cuvette into photometer



9. Time is displayed, wait for measured value

BIL / BIL N (BIL 142)



Sampling of BIL N

Additionally required: Mini centrifuge, Microvette




1. After pricking with the lancet, take about 60 μ L of blood (approx. 1 drop) from the heel with the microvette

Note: Close the microvette carefully before placing it into the mini centrifuge



2. Insert microvette into mini centrifuge
Centrifugate for 3-5 minutes
Note: Ensure an even loading inside the mini centrifuge

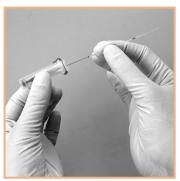

3. Take 20 μL plasma from the microvette

Continue with Fig. 1 on the previous page

Mini centrifuge

Art. No. DZ 002

Microvette

Art. No. LH 031 (100 pieces per package)

DIA QS / BIL QS



Quality assurance

Photometer testing with control caps







1. 20 control caps with lyophilized control serum

DIA QS: Glucose, Lactate, Triglycerides, Cholesterol, HDL-Cholesterol, Uric acid, Bilirubin (adults)


BIL QS: Bilirubin (adults), Bilirubin N (newborns)




2. Screw the control cap onto a cuvette of the test you want to check

Mix wel

Now the cuvette contains a sample with a known concentration

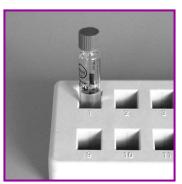

3. Leave the cuvette for 1 minute

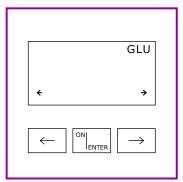

4. Switch photometer on with ON/ENTER key

Wait for device check, confirm with ON/ENTER

Select the required test, confirm with ON/ENTER



5. Zero point adjustment: Insert cuvette with sample (Fig. 3) into photometer, zero point is stored in memory Remove cuvette after signal tone



6. Replace control cap with starter cap



7. Turn cuvette upside down several times



8. First press ON/ENTER, then insert cuvette into photometer Wait for measured value Compare the measured value with the target value on the package insert